


μ FR[®] _ AngioPlus Core

Angio-based FFR and PCI Planning System



Pulse Medical

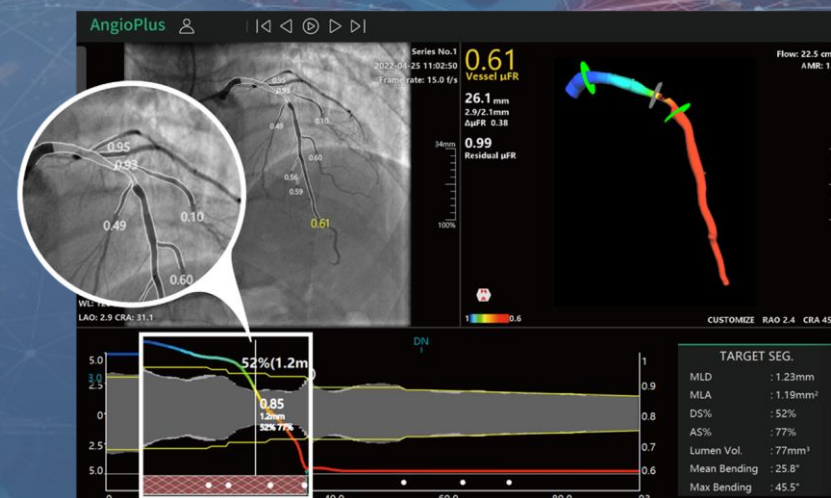

博动 PULSE MEDICAL

Clinical Using Scenario



France



Brazil

Clinical Using Scenario



Computer configuration



μFR (AngioPlus Core) Workstation configuration requirements

AngioPlus Core	Operating System	Microsoft Windows 10/11, 64bit
	CPU	Intel Core i7 Series 10th Gen or higher, or i5 Series 11th Gen or higher
	RAM	32 GB or higher
	Storage	At least 100 GB of available storage
	Monitor resolution	2560 * 1600, 2560 * 1440, 2880 * 1800, 1920 * 1080, 1920 * 1200, 3840 * 2160, 1600 * 1200
	Network card	100 Mbps or higher
	GPU (*optional)	<i>NVIDIA RTX 20 Series and 20 SUPER, NVIDIA RTX 30 series or NVIDIA RTX A4000, and NVIDIA RTX 4060/4070 (The GPU Memory should be at least 6GB).</i>
	GPU Driver (*optional)	<i>GeForce 471.11 Driver or later</i>

Product Package



Product Sticker

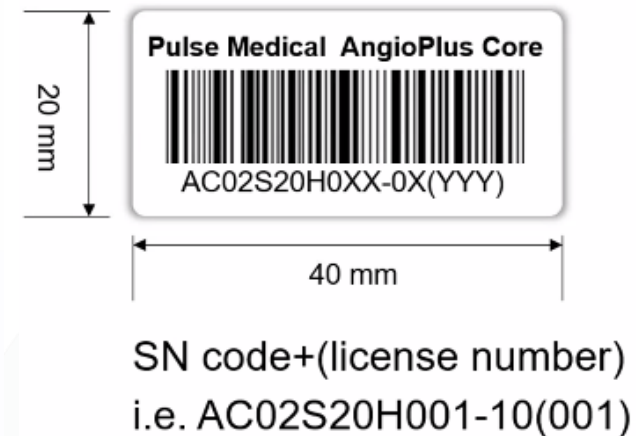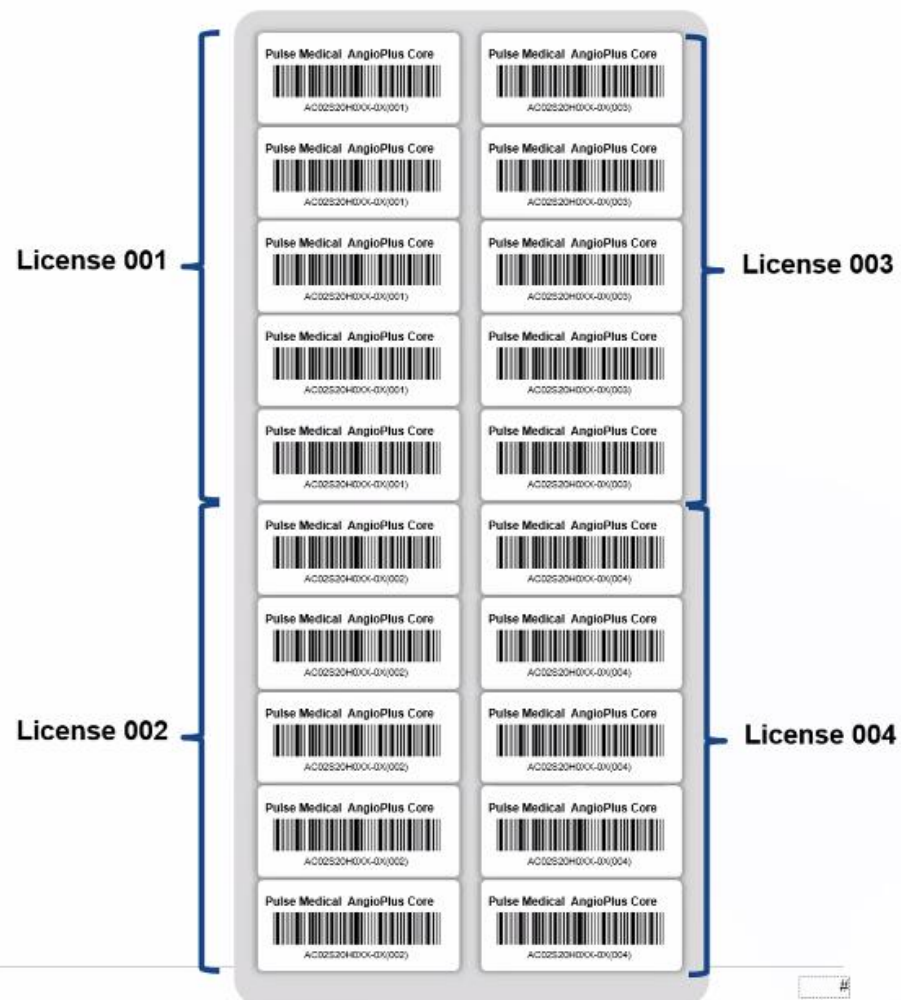

SN code+(license number)
i.e. AC02S20H001-10(001)

Analysis SOP

SOP of Coronary Angiographic Image Acquisition

1. Acquiring the angiographic image using recommended projection

Acquire **at least 1 angiographic image** runs according to the **recommended angiographic projections**. Adjust angiographic projections if the recommended projection has a suboptimal quality or vessel overlaps at the lesion segment.

2. Uniformly injection, fill the entire vessel

keep catheter to have better **co-axiality** with respect to vessel, and maintain the **uniformly** contrast injection speed (~4 ml/s) until the distal vessel is **filled** with contrast and keep it to include the **ED phase**.

3. Nitroglycerin required

Inject **nitroglycerin** through the imaging/guiding catheter to eliminate arterial spasm prior to angiography.

4. Angiographic frame rate ≥ 15 frame/s, start recording before contrast injection, ≥ 1 cardiac cycle

Set angiographic imaging frame rate **≥ 15 frame/sec**. Then prefill the catheter with contrast medium. After complete wash out of the contrast medium, Start recording angiographic images at least **1 second prior to contrast injection**. The angiogram should include **≥ 1 cardiac cycle**, preferably 2 cycles.


5. Table moving

Don't move the table, especially the height.



Analysis process

Image selection: Recommended projection, clear angiography, less overlapping, less foreshortening

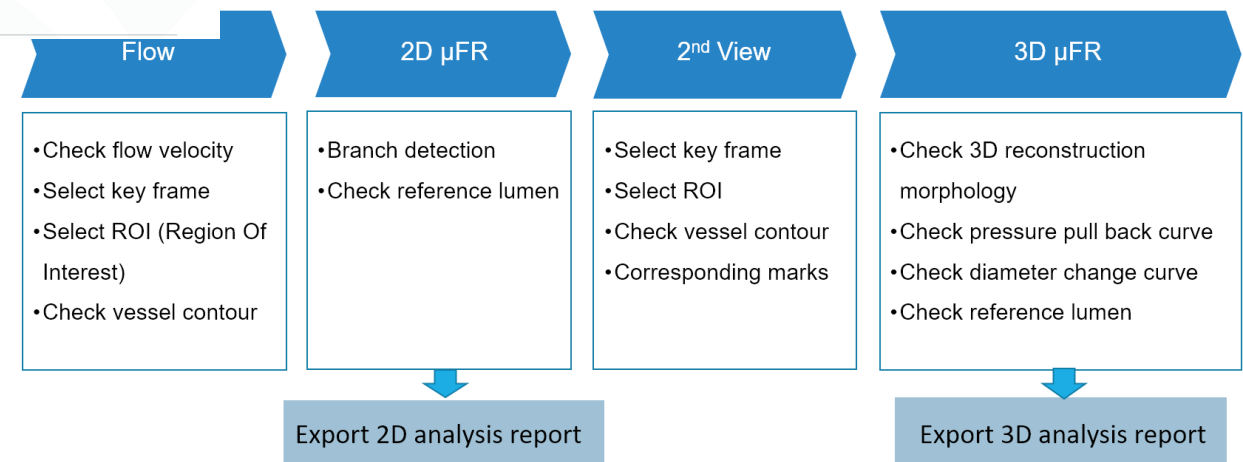

Data Connection

